NLRP1 Rabbit pAb

Catalog No.: A16212 3 Publications



Basic Information

Observed MW 160kDa

Calculated MW 166kDa

Category Polyclonal Antibody

Applications WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity Human, Mouse, Rat

Recommended Dilutions

Background

This gene encodes a member of the Ced-4 family of apoptosis proteins. Ced-family members contain a caspase recruitment domain (CARD) and are known to be key mediators of programmed cell death. The encoded protein contains a distinct N-terminal pyrin-like motif, which is possibly involved in protein-protein interactions. This protein interacts strongly with caspase 2 and weakly with caspase 9. Overexpression of this gene was demonstrated to induce apoptosis in cells. Multiple alternatively spliced transcript variants encoding distinct isoforms have been found for this gene, but the biological validity of some variants has not been determined.


Immunogen Information

WB	1:500 - 1:1000	Gene ID	Sw	iss Prot
	1 50 1 000	22861	Q9	2000
IHC-P	1:50 - 1:200			
IF/ICC	1:50 - 1:200	Immunogen Recombinant fusion protein containing a sequence corresponding to amino acids 1370-1473 of human NLRP1 (NP_127497.1).		
		Synonyms NAC; JRRP; MSPC; AIADK; CARD7; CIDED; NALP1; SLEV1; DEFCAP; PP1044; VAMAS1; CLR17.1; DEFCAP-L/S; NLRP1		
Contact	Product Information			
€	www.abclonal.com	Source	Isotype	Purification
		Rabbit	IqG	Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.


Validation Data

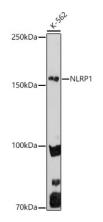


Western blot analysis of lysates from K-562 cells, using NLRP1 Rabbit pAb (A16212) at 1:1000 dilution.

Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021). Exposure time: 180s.

Immunohistochemistry analysis of paraffin-embedded human colon using NLRP1 Rabbit pAb (A16212) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

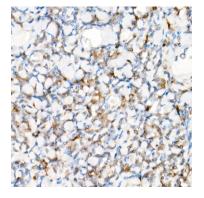


Immunohistochemistry analysis of paraffin-embedded human tonsil using NLRP1 Rabbit pAb (A16212) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

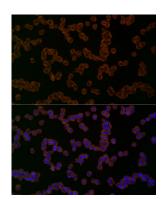

Immunofluorescence analysis of K-562 cells using NLRP1 Rabbit pAb (A16212) at dilution of 1:50 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.



Immunofluorescence analysis of PC-12 cells using NLRP1 Rabbit pAb (A16212) at dilution of 1:50 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

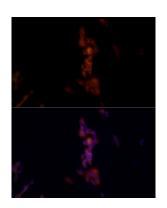

Immunofluorescence analysis of SH-SY5Y cells using NLRP1 Rabbit pAb (A16212) at dilution of 1:50 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.