Leader in Biomolecular Solutions for Life Science

PPAP2B Rabbit pAb

Catalog No.: A15743 1 Publications



Basic Information

Observed MW

35kDa

Calculated MW

35kDa

Category

Polyclonal Antibody

Applications

WB, ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

The protein encoded by this gene is a member of the phosphatidic acid phosphatase (PAP) family. PAPs convert phosphatidic acid to diacylglycerol, and function in de novo synthesis of glycerolipids as well as in receptor-activated signal transduction mediated by phospholipase D. This protein is a membrane glycoprotein localized at the cell plasma membrane. It has been shown to actively hydrolyze extracellular lysophosphatidic acid and short-chain phosphatidic acid. The expression of this gene is found to be enhanced by epidermal growth factor in Hela cells.

Recommended Dilutions

WB

1:500 - 1:1000

Immunogen Information

Gene ID 8613 **Swiss Prot**

014495

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human PPAP2B (NP 003704.3).

Synonyms

LPP3; VCIP; Dri42; PAP2B; PPAP2B

Contact

•

www.abclonal.com

Product Information


Source Rabbit **Isotype** IgG **Purification**Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

Validation Data

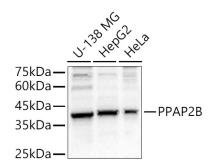

Western blot analysis of various lysates, using PPAP2B antibody (A15743) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.