

A0277

Leader in Biomolecular Solutions for Life Science



TNF- α Rabbit pAb

Catalog No.: A0277 **58 Publications**

Basic Information

Observed MW

Refer to figures

Calculated MW

26kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human,Mouse,Rat

Background

This gene encodes a multifunctional proinflammatory cytokine that belongs to the tumor necrosis factor (TNF) superfamily. This cytokine is mainly secreted by macrophages. It can bind to, and thus functions through its receptors TNFRSF1A/TNFR1 and TNFRSF1B/TNFR. This cytokine is involved in the regulation of a wide spectrum of biological processes including cell proliferation, differentiation, apoptosis, lipid metabolism, and coagulation. This cytokine has been implicated in a variety of diseases, including autoimmune diseases, insulin resistance, psoriasis, rheumatoid arthritis ankylosing spondylitis, tuberculosis, autosomal dominant polycystic kidney disease, and cancer. Mutations in this gene affect susceptibility to cerebral malaria, septic shock, and Alzheimer disease. Knockout studies in mice also suggested the neuroprotective function of this cytokine.

Recommended Dilutions

WB	1:500 - 1:2000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID

7124

Swiss Prot

P01375

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 77-233 of human TNF (NP_000585.2).

Synonyms

DIF; TNFA; TNFSF2; TNLG1F; TNF-alpha; TNF- α

Contact



www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.


Validation Data

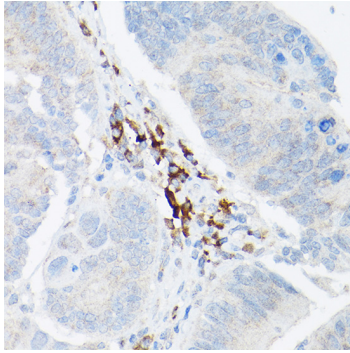


Immunohistochemistry analysis of paraffin-embedded human colon carcinoma using TNF- α Rabbit pAb (A0277) at dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

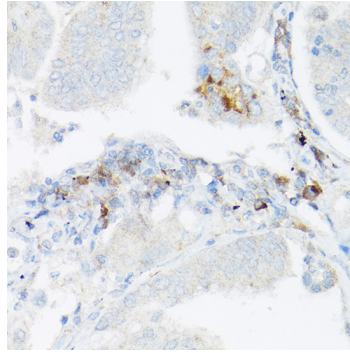

Immunohistochemistry analysis of paraffin-embedded human lung cancer using TNF- α Rabbit pAb (A0277) at dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

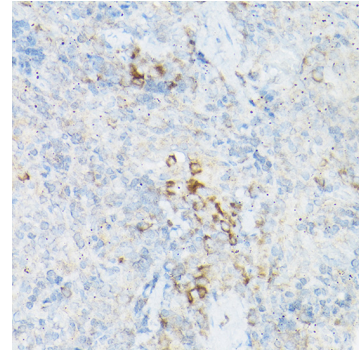

Immunohistochemistry analysis of paraffin-embedded rat spleen using TNF- α Rabbit pAb (A0277) at dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

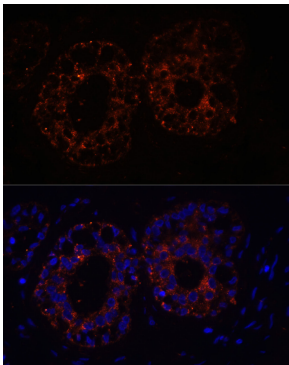

Immunofluorescence analysis of Human mammary cancer cells using TNF- α antibody (A0277) at dilution of 1:100. Blue: DAPI for nuclear staining.